



SPEEDY 4all Ergo

MOBILITY,
THAT ADAPTS TO YOU.

The **SPEEDY 4all Ergo** rigid frame wheelchair that is characterised by its fine frame geometry with ergonomic seat fold. With 9 frame parameters, it can be optimally adjusted to you. Many settings can also be made retroactively: The tipping point, seat heights or seat angle can be continuously adapted to your needs.




#Reliability

Sturdy and durable frame for maximum mobility and compatibility with auxiliary drives and adaptations.



#MoreMobility

Traction devices and adaptive bikes can be connected in just a few steps via adapters; optional wheelbase extension for extra traction.

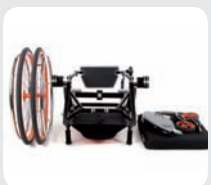

#Independence

Easy folding of the backrest and tool-free removal of assemblies enable independent transport.



#Flexibility

Integrated tipping point and seat height adjustment in 1 cm steps for optimal adaptation to individual needs.



PERFECTLY TAILORED TO YOU.

The SPEEDY 4all Ergo is a customised model according to your needs. You can choose from a wide range of equipment, including footrests, brakes, push handles, back and seat systems, clothing guards, armrests, as well as drive and caster wheels.


UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



The 7-way angle-adjustable and foldable backrest offers maximum flexibility for the ideal sitting position.



The footrest, which can be folded towards the back with a spring locking mechanism, is easy to handle and flexible in everyday life.

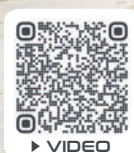

Manoeuvrable one-hand drive with double handrim for users with hemiplegia or restricted mobility on one side of the body.



Elegant integral parking brake with comfortable and safe operation.




TECHNICAL SPECIFICATIONS

Seat width:	30 - 52 cm (in 1 cm steps)
Seat depth:	36 - 48 cm (in 1 cm steps)
Back height:	20 - 48 cm (in 1 cm steps)
Lower leg length:	37 - 50 cm (in 1 cm steps)
Rear seat height:	38 - 47 cm (in 1 cm steps)
Front seat height:	44 - 52 cm (in 1 cm steps)
Tipping point:	3 - 11 cm (in 1 cm steps)
Total product weight:	Starting at 9.7 kg
Max. standard load weight:	Up to 120 kg
Max. load weight for Heavy Duty:	Up to 150 kg appropriate individual dimensions upon request

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

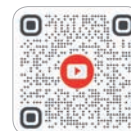

PRO ACTIV®

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 09/2024

Subject to technical modifications, printing errors & colour deviations.